




Z-7000 Series

STRAIGHT ZOOM STEREO MICROSCOPE

Z-7100 Straight Zoom Stereo Body

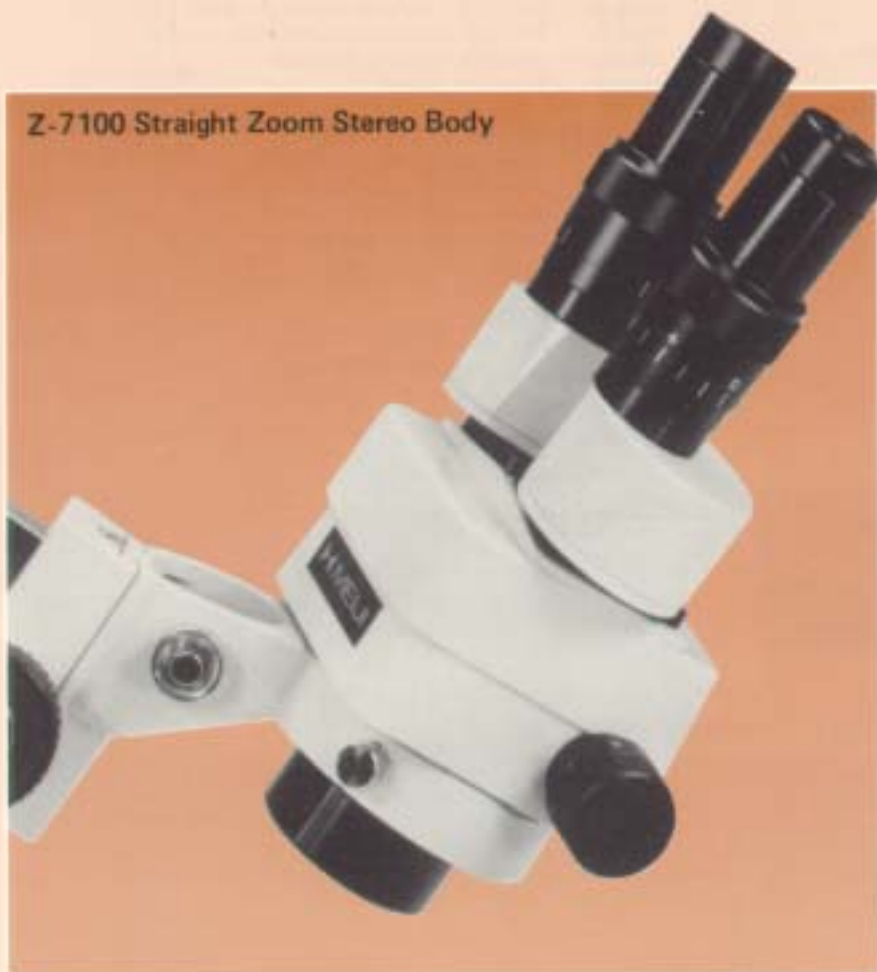

MEIJI TECHNO Z-7000 Series Straight Zoom Stereo Body forms eyepiece tubes and objective lenses straightforward in a line, focusing on the specimen directly through. Therefore, you can observe the object, actually facing it.

Another feature is in its extensible working distance which can be expanded as long as 330mm with the aid of an optional auxiliary lens.

With these features, MEIJI Z-7000 Series Straight Zoom Stereo Body can be an optimum microscope body to be installed on machines of various sorts in industry as well as the medical field where the long distance observation is often required.

MEIJI TECHNO Z-7000 Series Straight Zoom Stereo Body is mounted on machines easily by MEIJI F or FS Holder with a focusing block. It can surely expect a wider working area for efficiency.

F HOLDER

Designed to hold the zoom stereo body on machines or in special applications. It contains its own focusing slide block.



FS HOLDER

Designed to hold the zoom stereo body on machines or in special applications. Arm can incline through 180° and also pivot through 360°. It contains its own focusing slide block.



Supplementary lenses



ISSCO

Ph. (02) 9743 2588
Fax.(02) 9736 2835